

CareSens SmartLog Desktop: Adding a new profile (patient) using meter and USB cable


1. Plug the meter and cable in, turn the meter on (meter will say PC).
2. Click 'Download Data' on the left menu, and then the '**Download Data**' image.

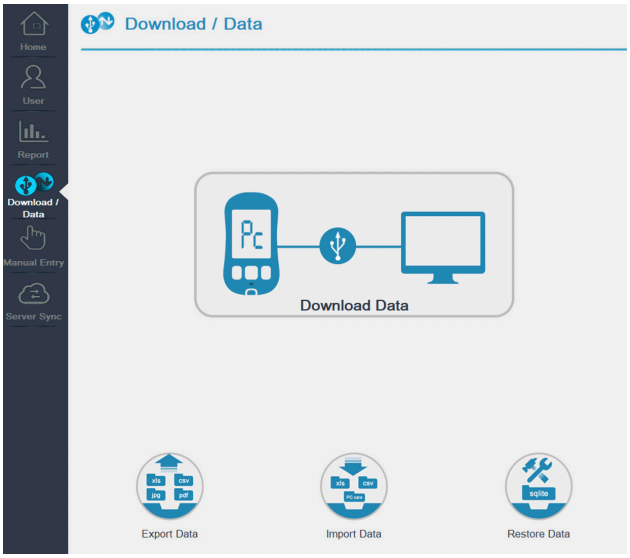


3. Press the disk/tick box (it goes green when you put the mouse over it).

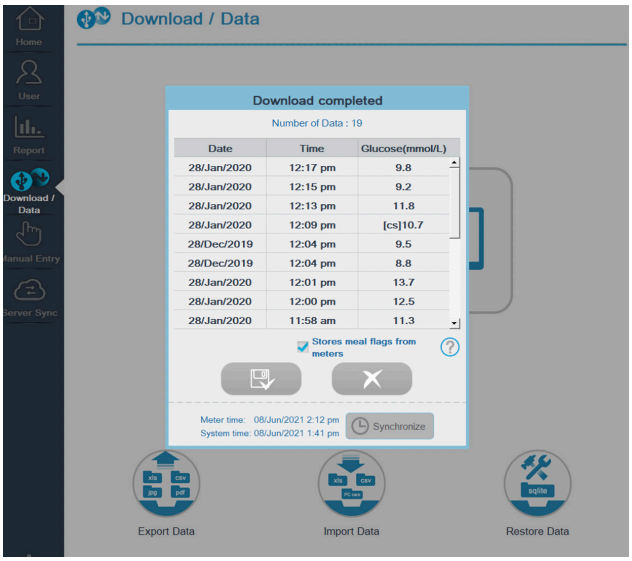

4. SmartLog detects and displays the meter serial number: for a new patient, leave the 'Name' dropdown box as '**Add user accoun**t', fill in their details, and press the tick box (bottom right, it goes green when you put the mouse over it).

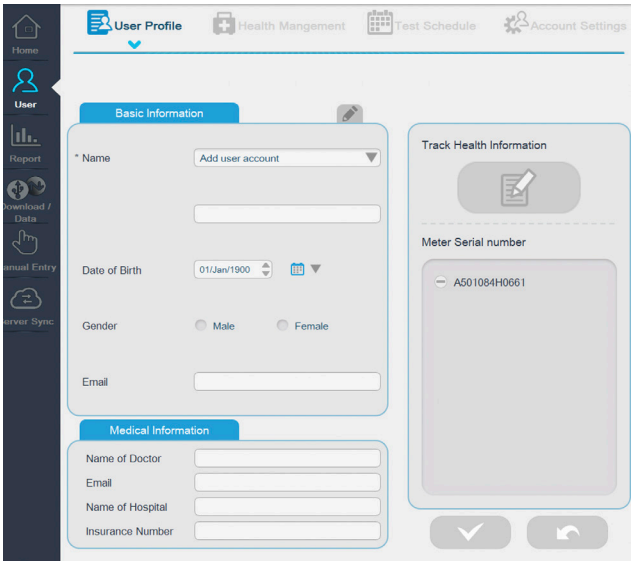


5. → add patient details and tick to save.

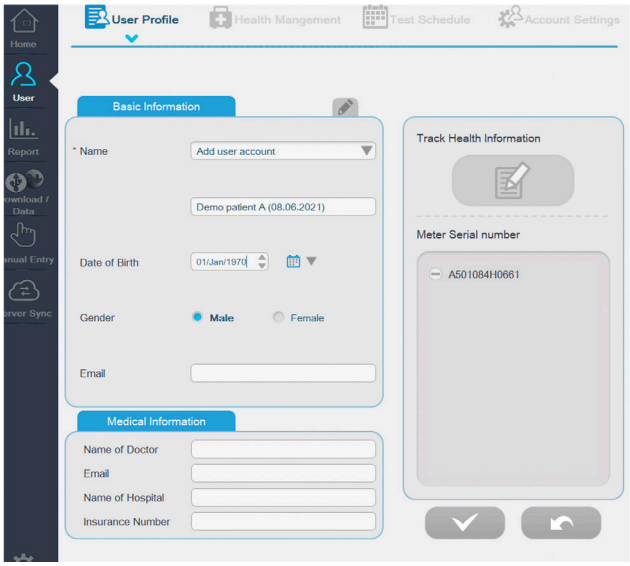


Pharmaco Diabetes


6. A message appears that new data has been successfully added, and they appear in the patient/profile list at the top of the screen.

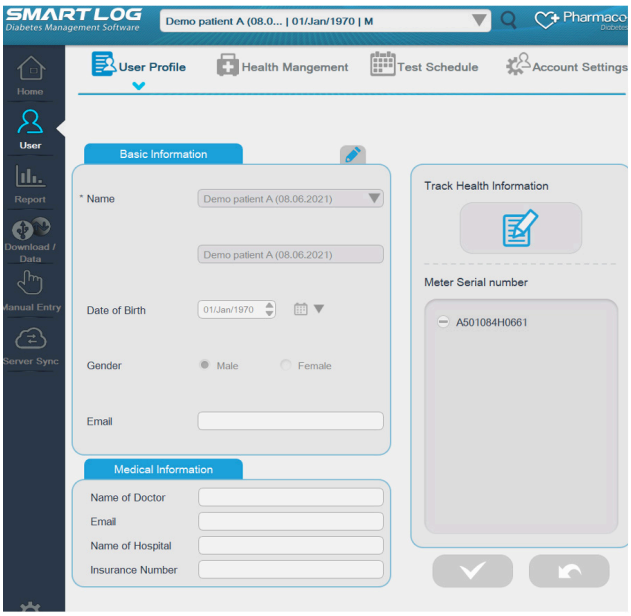


7. Repeat for the next new patient – plug their meter in, download data and set up their profile.

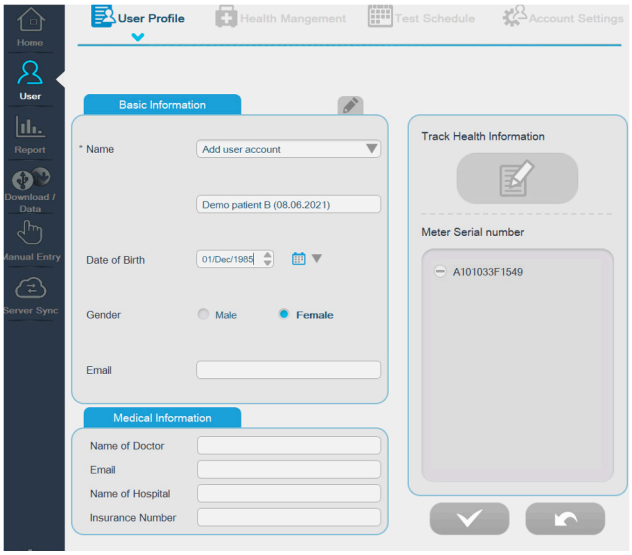


8. You can see the patients registered in either the drop-down box at the top of the screen, or by going into '**Account Settings**'.

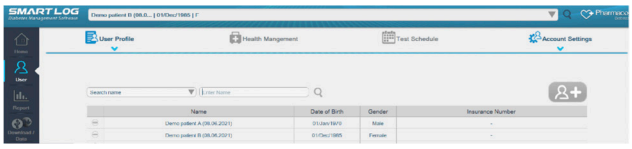


9. Next time the patient comes in (they already have a profile set up), plug in the meter, download and SmartLog will upload the data to their profile, as it recognises their meter serial number.
e.g., Patient A previously set up, and the first import had these results:

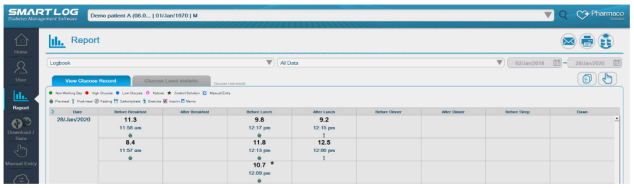

10. Upload the meter today, which will add data since the last upload: plug the meter in, download data as above – SmartLog recognises the meter serial number, and automatically brings up the patient profile registered to that serial number, and imports the meter readings to Patient A. A message will display that new data has been successfully added.

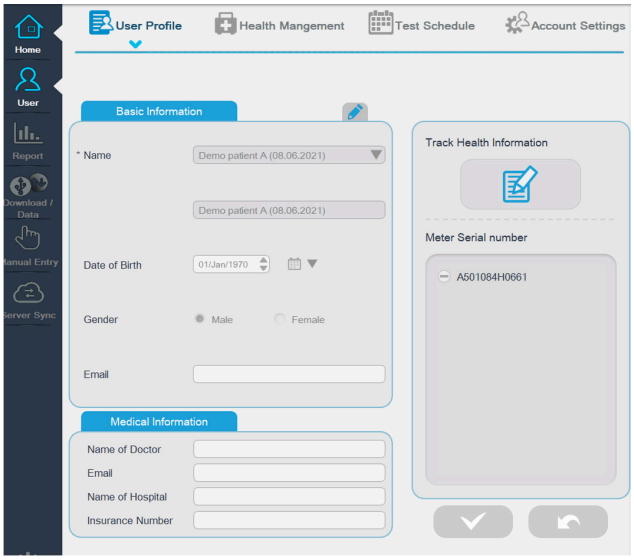

11. Go into the report – these more recent readings have been added.

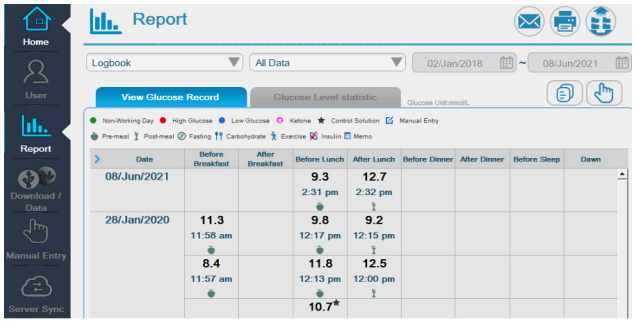

12. The manual can also be downloaded from SmartLog.

